
Smartline

▶ Column Oven 4050 Manual/Handbuch

V7335A / V7335

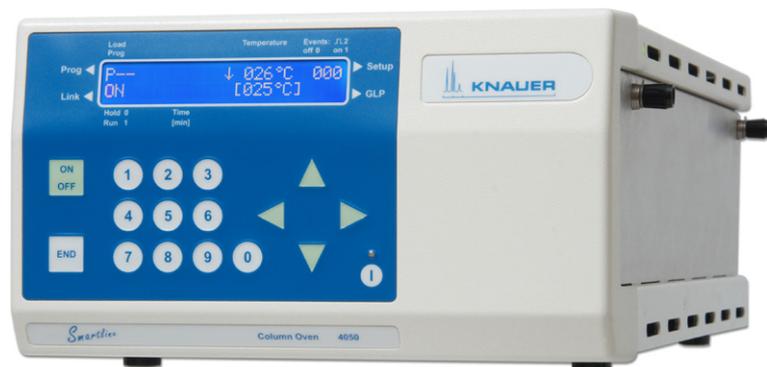Contents

Note: For your own safety, carefully read the manual and always observe the warnings and safety information on the device and in the manual!

Contents	1
Using this Manual	5
Conventions in this manual	5
SOP's in this manual	5
Intended use	6
Device types	6
Laboratory use	6
Where is it prohibited to use the device or system?	6
Safety	6
Laboratory regulations	6
Solvents	7
PEEK connections	8
Protective measures	8
Power supply and mains connection	8
Target group	8
Operating the device or device system	8
To what should the user pay particular attention?	9
What expertise should users have to safely operate a HPLC device or device system?	9
Symbols and labels	9
Installation	10
Packaging and transport	10
Fastening material and shipping boxes	10
Protective film on the display	10
Scope of supply	10
Checking the scope of delivery	11
Space requirements	11
Installation site	11
General Description	12
Operating the Column Oven 4050	13
Unpacking	13
Function Elements of the Column Oven 4050	13
Front Door of the Column Oven 4050	13
Function of Foil Keys	14
Front elements of the Smartline Column Oven 4050	15
Side view of the Smartline Column Oven 4050	15
Rearview of the Smartline Column Oven 4050	16
Power Supply, ON/OFF	16
Internal Software Structure	17
Setup Screens	18
Prog Screens	19
Program Editing Screens	19
Deleting Programs	20
Link Editing Screens	21
GLP Menu	22
ERROR Screens	22
Position of the Column Oven 4050 in Smartline Systems	23
Installation of the Columns	23
Access to the Column Compartment	23
Column and Capillary Holding	24
Direct Control of the Column Oven 4050	27
Software Control of the Column Oven 4050	27
RS-232 Serial Interface	28

LAN – Ethernet Interface	28
Analog Interface	29
Using the remote control socket	29
Connections to the remote control socket	29
Configuring Plug Strips	30
Installation Qualification (IQ)	31
Operation Qualification (OQ)	31
Test intervals	31
Execution	31
Maintenance and care	31
Technical support	31
Maintenance contract	31
What maintenance tasks can users perform on the device?	32
Cleaning and caring the device	32
Environmental protection	32
Disposal	32
Decontamination	32
Storage	32
Ambient storage conditions for the device	32
Technical Data	33
Delivery program	34
Accessories	34
Spare parts	35
Legal information	36
Warranty conditions	36
Manufacturer	36
Declaration of conformity	37
Index	38

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Using this Manual

This manual refers to the Smartline Column Oven 4050.

Conventions in this manual

Arrows like this: \longleftrightarrow , \rightarrow , used in block diagrams, indicate that the user is asked to press the corresponding arrow keys. The operation of arrow keys is defined as follows:

Cursor up: ▲

Cursor left: ◀

Cursor right: ▶

Cursor down: ▼

Important Hints are marked by the marginal hand symbol.

Special warnings are indicated by the marginal warning sign and printed in bold letters.

The marginal lamp symbol indicates helpful advice.

SOP's in this manual

The **Standard Operating Procedures (SOP)** provided with this manual offer a convenient way of structuring complex tasks in the operation of your Column Oven 4050. They include step-by-step instructions leading the user through all routine tasks during operation. They can be used for documentation purposes and be copied, applied signed, and filed in order to document the performance of the instrument.

Please operate the instrument and all accessories according to instructions and SOP's in this manual. This ensures accurate results and prolongs the longevity of your equipment.

SOP 1	Opening and closing the column compartment	23
SOP 2	Installation of the capillaries	24
SOP 3	Installation of the column	24
SOP 4	Mounting of equipped column holdings	25
SOP 5	Connecting Capillaries	26
SOP 6	Mounting the outlet capillary into the cooler	26
SOP 7	Mounting Plug Strips	30

Please Note the Smartline Column Oven 4050 can reach a maximum temperature of 145°C. Normal silica gel containing HPLC columns are not allowed to be operated above 90°C. Please set the allowed maximum temperature for your columns in use. Also observe the temperature durability of all the other components that are used.

Intended use

High performance liquid chromatography (HPLC) is a method for separating substance mixtures, determining substances and measuring their concentration.

The device or system is suitable for high-pressure liquid chromatography. It is suitable for laboratory use, for analyzing substance mixtures that can be dissolved in a solvent or solvent mixture.

Device types

- Smartline Column Oven 4050

Laboratory use

- Biochemistry analyses
- Chiral analyses
- Food analyses
- Pharmaceutical analyses
- Environmental analyses

Where is it prohibited to use the device or system?

Risk of explosion! Never use the device in potentially explosive atmospheres without appropriate protective equipment and approval by a notified body! Inform the technical support department of the manufacturer!

Safety

Laboratory regulations

Observe national and international regulations pertaining to laboratory work!

- Good Laboratory Practice (GLP) of the American Food & Drug Administration
- For development of methods and validation of devices: Protocol for the Adoption of Analytical Methods in the Clinical Chemistry Laboratory, American Journal of Medical Technology, 44, 1, pages 30–37 (1978)
- Accident prevention regulations published by the accident insurance companies for laboratory work

Solvents

Even small quantities of other substances, such as additives, modifiers, or salts can influence the durability of the materials. The list of selected solvents was compiled based on research in the pertinent literature and is only a recommendation by the manufacturer. In the event of doubt, contact the technical support of the manufacturer.

Solvents suitable for use in HPLC:

- Acetone
- Acetonitrile
- Benzene
- Chloroform
- Acetic acid (10–50%), at 25 °C (77 °F)
- Ethyl acetate
- Ethanol
- Hexane/heptane
- Isopropanol
- Methanol
- Phosphoric acid
- Toluol
- Water

To avoid damage from leaks, always place solvent bottles in a solvent tray on the device.

Organic solvents are toxic above a certain concentration. ensure that work areas are always well-ventilated! Wear protective gloves and safety glasses when working on the device!

Organic solvents are highly flammable. Since capillaries can detach from their screw fittings and allow solvent to escape, it is prohibited to have any open flames near the analytical system!

Only use solvents that have a self-ignition point higher than 150 °C under normal ambient conditions!

The following solvents can attack the components of the device and are therefore not suitable:

- Mineral and organic acids (except in buffer solutions)
- Bases (except in buffer solutions)
- Liquids containing particles

The following solvents are suitable to only a limited extent for use in the device:

- Methylene chloride
- Tetrahydrofuran (THF)
- Dimethyl sulfoxide (DMSO)
- Slightly volatile solvents
- Fluorinated hydrocarbons

Regularly check for leaks and clogged capillaries – test back pressure without column!

PEEK connections

- Disposable PEEK fittings
- One-piece disposable polyetheretherketone fittings (PEEK) for easier installation of flexible, thin capillaries
- Tightening torque of the PEEK screw: Tightened by hand (approx. 0.5 Nm)

Protective measures

- You are only permitted to perform the maintenance tasks described in this manual.
- All other maintenance tasks are to be performed exclusively by the manufacturer or a company authorized by the manufacturer.

Without exception, the following applies to all maintenance tasks that can be performed by the user:

- Switch off the device and pull the power plug!
- Never open the device! High voltage poses a life-threatening risk!

Power supply and mains connection

The device is intended for use with AC power networks of 100–240 V. The supplied power cable is to be used to connect the device to the mains supply.

To disconnect the device from the mains power, disconnect the power plug on the rear of the device.

Target group

Operating the device or device system

The device can be operated by means of the chromatography software at the workstation or the membrane keyboard at the device.

To what should the user pay particular attention?

To make your HPLC separations as efficient as possible, pay close attention to the following:

- Once they have been used, never re-use capillaries in other areas of the HPLC system.
- Only use a given PEEK fitting for one specific port and never re-use it for other ports. Always install new PEEK fittings on each separate port.
- Use special columns and follow the manufacturer's instructions on caring for the columns!
- Regularly check for clogged capillaries – test back pressure without column!
- Use ultra-pure, filtered solvents – Gradient grade – for the HPLC.
- Filtration of substances to be analyzed.
- Use of inline filters.

Only allow the technical support department of the manufacturer or a company authorized by the manufacturer to open the devices for maintenance and repair work.

What expertise should users have to safely operate a HPLC device or device system?

- Completed degree as chemical laboratory technician or comparable vocational training
- Fundamental knowledge of liquid chromatography
- Participation in an installation of the system performed by the manufacturer or a company authorized by the manufacturer, or suitable training on the system and chromatography software
- Basic knowledge of Microsoft Windows®
- Knowledge regarding substances that are suitable only to a limited extent for use in liquid chromatography

Symbols and labels

Explanations of symbols and labels on the device or system:

Hazard symbol indicating microelectronic devices that can be damaged by electrostatic discharge when touched.

CE-marking for devices that comply with the pertinent EU directives (Conformite Europeenne) and come with a declaration of conformity from the manufacturer.

For your own safety, read the manual and always observe the warnings and safety information on the device and in the manual!

Marking for devices that comply with the Canadian requirements for laboratory equipment: CAN/CSA-C22.2 No. 61010-1, second edition, including Amendment 1, or a later version.

Installation

Packaging and transport

At the factory, the device was carefully packed for safe transport.

Check the device for signs of damage that occurred during transport. If the shipment is incomplete or damaged, inform the manufacturer within three workdays. Also inform the freight carrier about transport damage.

Fastening material and shipping boxes

The device is held in place and protected by foam inserts at the top and bottom. Please keep the transport box and the foam inserts.

1. Remove the foam insert on the top of the device.
2. Grip the device at its sides, near the front panel, and lift it out of the packaging.

Protective film on the display

During transport, a protective film prevents scratches to the display of the device.

Remove the protective film from the display!

Scope of supply

Column Oven 4050

- Manual
- 2 x half shell holdings for columns of 9 to 12 mm OD 100 mm
- 4 x half shell holdings for columns up to 8 mm OD 100 mm
- Fixing plate 100 mm
- Fixing plate 200 mm
- 2 x fixing plate 40 mm
- 2 x half shell holdings for columns of 9 to 12 mm OD 40 mm
- 2 x half shell holdings for columns up to 8 mm OD 40 mm
- Set WAGO Micro connect. strips for Smartline 8pol.
- 10-polar flat band cable
- 16 x ISO 4762 M3 x 10 A2 cylinder screw
- 8 x ISO 4762 M3 x 20 A2 cylinder screw
- Hexagonal key no. 2,5
- Hexagonal key No. 3
- Double open-end wrenches, 1/4" and 5/16"
- Power Cable type VIIG C5
- Null Modem cable (complete)
- Patch cable CAT 5e

Use original parts and accessories made by the manufacturer or a company authorized by the manufacturer.

Checking the scope of delivery

1. Check whether the device and accessories are complete.
2. If anything is missing, inform the technical support department of the manufacturer.

Hotline of the technical support department of the manufacturer:

Languages: German and English

Available by telephone: 8 am to 5 pm (CET)

Phone: +49-(0)30-809727-0

Fax: +49-(0)30-8015010

email: info@knauer.net (manufacturer)

Space requirements

Side clearance to other devices:

- If there is a device on one side, min. clearance of 5 cm.
- If there are devices on both sides, min. clearance of 10 cm.
- At least 30 cm gap to the fan on the rear of the device.

Make sure that the power plug on the rear of the device is always accessible, so that the device can be disconnected from the power supply.

Installation site

- Air humidity: Below 90% (non-condensing)
- Temperature range: 4–40 °C; 39.2–104 °F
- Sunlight: When setting up the device at the installation location, make sure that it is protected against direct sunlight.

General Description

Fig. 1 *Smartline Column Oven 4050*

The **Smartline Column Oven 4050** allows for maximum stability of your analytical measurements by preventing temperature changes from having an effect on your HPLC columns. A stable column temperature not only improves reproducibility, but considerably improves the resolution of many HPLC methods as well. Thanks to its compact design, this column thermostat fits perfectly into every HPLC system of the Smartline series yet can also be conveniently integrated in other systems. All functions of the oven can be easily controlled with the HPLC software packages ChromGate[®] or ClarityChrom[®] as well as directly from the front panel.

The **Smartline Column Oven 4050** works very effectively as a contact oven while also including thermostating of the mobile phase at the column inlet. This ensures an even temperature distribution over the entire column. Besides offering accurate temperature control for standard applications, high temperature HPLC and GPC applications up to 140 °C are also supported. Precise temperature values are made possible over the full operational range by using optimized components for measurement and control.

The **Smartline Column Oven 4050** features heating and active Peltier cooling of the column compartment and is suitable for applications above and below ambient temperature. With its temperature range from 5 °C up to 140 °C, the majority of HPLC applications can be realized. The oven also features an active cooling of the column eluate. This cooling device provides for optimal results in the subsequent detection through minimized noise and drift values.

The **Smartline Column Oven 4050** can be equipped with two columns up to 12 mm OD and a length of 300 mm or up to 4 columns of the same OD and a length of 125 mm. A mixed mounting of one long column and one or two short ones is also possible.

Operating the Column Oven 4050

Unpacking

All KNAUER instruments are always packed carefully and safely for transportation. After unpacking, please check the device and accessories thoroughly for any damage that may have occurred during transport. If necessary, put forward any claim for damages to the carrier.

Use the list "Scope of supply" on page 10 and check if the Column Oven 4050 is complete. Please contact your local dealer or the KNAUER headquarters (sales department) in Berlin, Germany if you are missing something or if you need support.

Remove the transport protection from the display.

Function Elements of the Column Oven 4050

Front Door of the Column Oven 4050

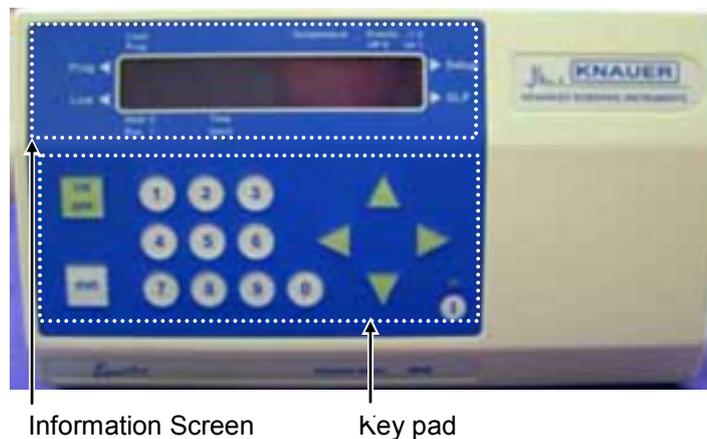
Fig. 2 Front view of the Column Oven 4050

Fig. 2 shows the front panel of the Column Oven 4050. The information screen and the key pad are arranged on the housing door. Capillaries are behind the door to prevent damaging.

Around the screen the links are placed to reach the menus Prog◀, Link◀, ▶Setup and ▶GLP.

Load and Prog, Temperature, Event off 0, on 1 und \square 2, Hold 0 und Run 1 as well as Time[min] are to describe the screen.

Function of Foil Keys

The keypad in Fig. 2 "Front view of the Column Oven 4050" consists of four control keys (arrow keys), two square keys (ON/OFF and END), and the standby key.

ON OFF - Key

This key is to switch ON/OFF the heating/thermostating. In all other menus pressing the ON/OFF key will result in an immediate return to the monitor screen.

END - Key


This key ends or interrupts the actual function and returns the monitor screen.

Numeric keys

The numeric keys are to enter data and to control program sequences.

Arrow Keys

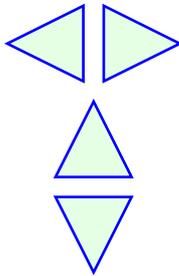
The arrow keys in light green are used for moving the cursor on the display and to confirm settings.

Use the arrow keys ► "right" or ◀ "left" to move from one menu parameter to the other as well as to confirm entered or selected values. Pressing the keys for a prolonged time changes the monitored menu.

Use the arrow keys ▲ "up" or ▼ "down" to change the selected parameter or its available options.

Standby Key

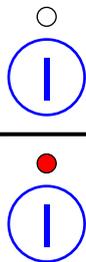Pressing the standby key for more than two seconds switches the instrument off (only the standby electronics will still be working). In standby mode the red LED integrated into the standby key is on. To switch the instrument on again, the standby key has to be pressed for more than one second. Then the instrument powers up again and the red LED is off.

Front elements of the Smartline Column Oven 4050

Behind the front door the separate cooling of the post column capillary is accessible. It can be used for 1/16" capillaries (slots for thermo transfer in the cover part) and for 0.5 mm capillaries (slots for thermo transfer in the inbuilt part).

Fig. 3 *Post column cooler*

To prevent the condensation of air humidity on the cooler it should be shielded with the isolation cap while the cooler is in use (Fig. 3).

Fig. 4 *Post column cooler with isolation cap*

Side view of the Smartline Column Oven 4050

The Smartline Column Oven 4050 can be equipped with two columns up to 12 mm OD and a length of 300 mm or up to 4 columns of the same OD and a length of 125 mm. A mixed mounting of one long column and one or two short ones is also possible. For shorter columns are corresponding column carriers available.

Fig. 5 shows the side view of the column oven 4050. The LED display at the housing shows the temperature state:

- Green LED: the entered "door" temperature is not exceeded.
- Red LED: the entered "door" temperature is exceeded.

Risk of burning!

The "door" temperature serves the safety settings. The default value is 60°C.

Side panel elements

- 1 LED display (red/green)

Fig. 5 Opened Oven compartment with a 125 mm column

Rearview of the Smartline Column Oven 4050

Rear panel elements

- | | | | |
|---|------------------------|---|---------------------------|
| 1 | ON/OFF Switch | 4 | RS232 Interface connector |
| 2 | Power connector | 5 | LAN connector |
| 3 | Events connector strip | 6 | Fan |

Fig. 6 Rearview of the Smartline Column Oven 4050

Power Supply, ON/OFF

The Column Oven 4050 is equipped with a modern switching power supply. This power supply works in a range of 90 - 260 V AC and 47-63 Hz. The instrument can be switched off and on again from the main switch (on the back of the instrument) or the standby key.

Please note that the instrument is not completely switched off when in standby-mode. Only with the main switch can the instrument be completely switched off.

Connect the Column Oven to the power supply cable and switch on the instrument. The power switch is on the rear panel. After switching ON, the display presents briefly information concerning the instrument and the firmware version:

**HPLC – COLUMN OVEN
V01.00**

After successfully completing the initial checks the oven is ready for use. In the display the main menu is shown:

**SYSTEMTEST
OK**

This is followed by the monitor screen on the display. The same situation is given if the oven is switched on from the standby status.

**P— ↑ 025°C 000
OFF [030°C]**

Above the ON/OFF status the active/loaded program, the actual temperature, and the status of external events are shown. The up or down arrow prior to the temperature value indicates if the actual temperature is below or above the set value, shown in brackets.

According to the setting (see setup screens below) the temperature is given in °C or °F. If the oven is within the temperature equilibrium the screen is simplified:

**P— 025°C 000
ON**

The set temperature will automatically be shown as the actual measured temperature is outside the range defined in the firmware and/or the cursor is moved to the temperature field.

If the cursor is moved to the field for the set temperature the value can be set from 15 °C up to 140 °C can be set using either ▲ or ▼ keys or directly the numeric keys.

With the cursor on the events field these can be defined with the 0, 1, or 2 key or using again the ▲ or ▼ keys. The meanings are 0 for OFF, 1 for ON and 2 for pulse ▭.

Following operating states are possible (see Direct Control of the Column Oven 4050 on page 27):

OFF / ON / RUN / HOLD / WAIT / REMOTE (⊗ ON) / REMOTE (⊗ OFF)

Internal Software Structure

The software is divided in various menus, each of which allows particular settings and operational modes.

You reach the single menus by positioning the cursor on the rhombus field ♦ followed by pressing the ▲ or ▼ keys. As indicated in Fig. 7, the menu will be called up in an endless loop.

Setup Screens

Move the cursor by shortly pressing the ► or ◀ key to the upper screen row. From here a prolonged pressing of the ► key will open the setup menu loop.

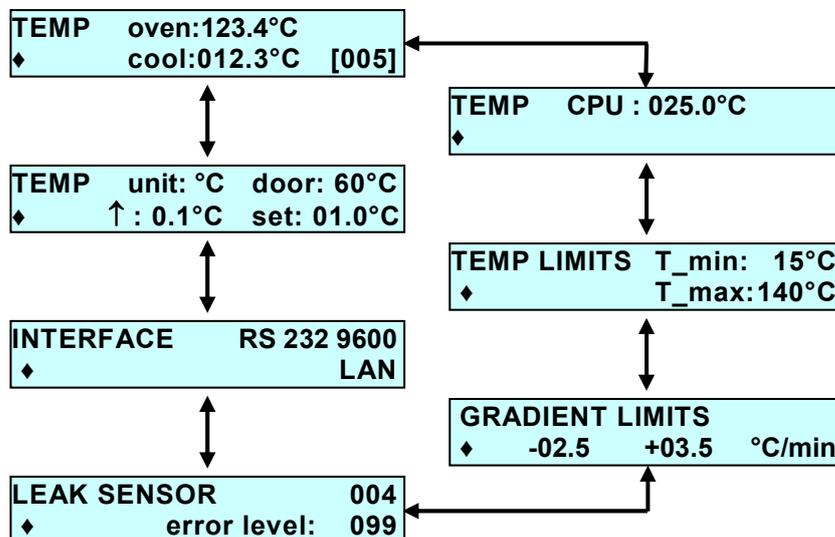
Fig. 7 Menu Sequence of the Column Oven 4050

The temperature monitor displays the actual temperatures of the oven and the capillary cooling. The set point can be entered within the range of 5 °C up to 140 °C.

The next screen displays the actual CPU temperature.

The temperature limits can be set within the instrument limits (5 °C through 140 °C) to deviating values. These settings then are limiting the possible set points. The minimum temperature of the oven is maximum 10 °C below room temperature.

In the Gradient limits menu, the minimum and maximum heat rates are set for the temperature program. The minimum heating or cooling rate is 0.1 °C per minute; the maximum heating or cooling rate is 10 °C per minute.

In the next screen of the temperature monitor the "door" temperature can be set. The temperature range is between +5°C and +100°C. The default value is +60°C.

The temperature will be set as follows: Move the cursor by shortly pressing the ► or ◀ key to the field "door: 60°C". The desired temperature will be set by either pressing the ▲ or ▼ keys or by pressing the number keys.

The next screen shows the actual value of the leak sensor and enables the editing of an error threshold. The threshold can be set from 1 (extreme sensitivity) through 100 (extreme low sensitivity, with the meaning OFF).

Setting the Interface:

The cursor blinks on the diamond key (◆). Move the cursor by shortly pressing the ► or ◀ key to the field RS 232 9600. Select LAN by pressing the ▲ or ▼ keys. Exit the menu by pressing the ► or ◀ keys.

The temperature unit can be switched through and back from °C to °F. This setting converts the all temperature values including the set points and limits.


On the ↑ field the interval for displaying the set points can be defined. Deviates the actual temperature more than this setting from the set point leading up or down arrows will indicate the direction of the deviation.

On the **set** field a second hysteresis can be set which influences the set point monitoring. Deviates the actual temperature more than this setting from the set point the main screen will display the set point in the second line in brackets.

Prog Screens

The editing menu are accessible by prolonged pushing the ◀ key. If the cursor is in this moment in the upper line, you will get the program editing menus. From the lower line you will be led to the link editing menus.

The prog menu consists of two screens and provides an overview on already existing and stored programs. Empty program places or numbers are indicated by "--". Here you can select the program number to be edited.

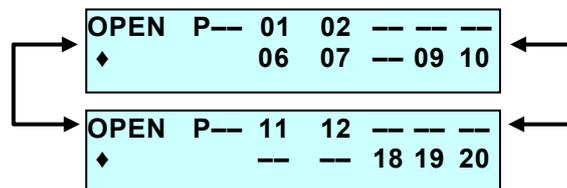There are 20 program sites available with a total of 100 program lines. Existing programs are indicated by their numbers and free places by "--".

Behind „Open P“ generally appears „--“. Enter the number for the program to be edited.

Running programs cannot be selected.

**ERROR Program is used
In a link !**

If a loaded but not running program is selected it becomes reloaded for editing.

The attempt to select a non existing program will lead to a corresponding error message.

**ERROR Program or link
Does not exist !**

Program Editing Screens

All program editing screens show in the lower left hand corner a rhomb (◆). Starting from this field the ▲ or ▼ keys lead to the following or preceding program line.

Programs to be edited are called up with their first line. The program number e.g. „P03“ is not free for editing. It is only information regarding the edited program.

New programs are called up with default values (25 °C and events 000) for the first line.

Parameter fields showing „-“ have a special meaning. These parameters do not influence the program flow!

Run Time Settings:

Run times are defined with three decimal places.

The first line is fixed to „000.0“. No entries can be made in this field of the first line. In any other line entering the time „000.0“ will delete this line.

Entering a new time or changing the time in any row will automatically bring the lines in the proper time order.

If the new time already exists in any other row this line will be overwritten.

Screens of Example Programs:

```
P01          000
◆ 000.0    025°C
```

```
P01          000
◆ 005.0    035°C
```

```
P01          000
◆ 010.0    035°C
```

```
P01          000
◆ 020.0    025°C
```

The screen for creating new program lines can always be found as the last row of the edited program:

```
P01          ** new time **
◆ 020.0
```

The time of the last line will be preset for the new line. After entering the new time the screen for line editing will appear:

```
P01          ---
◆ 025.0    ---°C
```

In case there is a faulty condition in the oven, the heater switches itself off.

Deleting Programs

If you enter „00“ in the prog menu a function for deleting programs is activated.

```
Delete program ?  P --
```

Enter the number of the program to be deleted. Entering „00“ once more will delete all existing programs. Otherwise only the program of the entered number will be deleted. In case there is no input, the display jumps back to the Prog screen (overview).

```
ERROR Programm or Link
does not exist!
```

If the input is wrong, an error message appears on the display. The error message can be exited with ▲ or ▼.

Just running programs or those which are integrated into a link cannot be deleted. The attempt brings up the corresponding error message:


```
ERROR Program is used
in a link
```

After an error message or successful deletion, there is a jump back to the Prog screen (overview).

Links can be deleted in the very same manner.

Link Screens

The link menu consists of two screens.

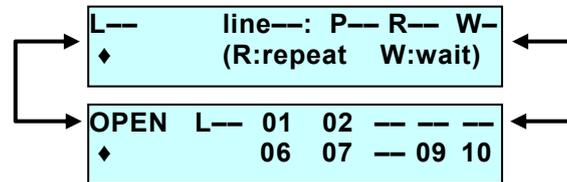The first one provides an status overview of the running link: Link no., actual line, actual program in the nth repeating, and the Wait status.

The second one provides an overview already existing and stored links. Empty link places or numbers are indicated by "--". Here you can select the link number to be edited.

Link Editing Screens

Structure and handling of the link editing screens are very similar to those of program editing. Instead of the times the row number is shown and used as order principle.

Screens of Example Links:

```

L03
◆ line01: P09 R15 W1
  
```

The screen for creating new link lines can always be found as the last row of the edited link. The number of the last line will be preset for the new line.

```

L01 ** new line **
◆ line 28
  
```

After entering the new number the screen for line editing will appear, adopting the settings of the previous link line:

```

L03
◆ line29: P17 R03 W1
  
```

New links are called up with default values (line 01, first available program, R01, W1) for the first line.

The line number counts continuously from the first line (01) to the last one. Entering a new number will move this line to the corresponding position. For instance entering 10 in the 5th line will move this line to position 10. Entering "00" will delete this line.

For the wait status two possibilities are available, W0 and W1.

W0 The first run of a program will wait for an external trigger signal. All following runs will proceed without delay. Any program start generates the trigger signal for the next run.

W1: All runs within the given link row will wait for the external trigger.

The deleting of links can be done in the same way as it was described for deleting programs, see page 20.

GLP Menu

The GLP menu is accessible by prolonged pushing the ► key while the cursor is placed in the lower line.

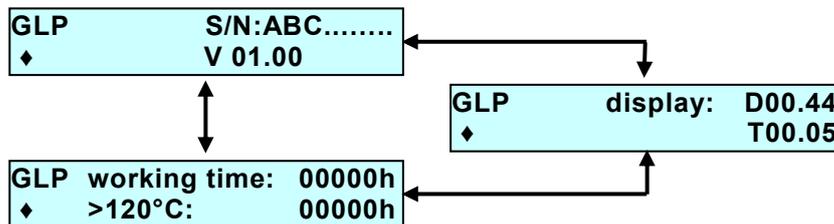
Fig. 8 GLP-Menu sequence

In the GLP menu the instrument's serial number and firmware version are indicated.

The display screen shows versions of internal firmware parts which are of interest only for service purposes.

The next screen monitors the operating time with the heating/cooling switched on as well as the operating time at high temperatures.

ERROR Screens

Besides the already described error messages related to the selection of links or programs (see page 19) the following ERROR messages can be occasionally displayed.

ERROR Leak detected !
◆ 234

Check the capillary connections!

ERROR Door not closed !

Close the oven door.

ERROR Programm or Link does not exist!

Wrong input during deleting of programs.

ERROR Temperature out of limits !

The temperature limits in the setup menu were changed. Calling up a previously written program with temperatures outside this new range brings up this message. Change according to the requirements either the temperature in the program or the set limits in the setup menu.

All ERROR messages will be deleted by pushing the ON/OFF or END key. In any case if a ERROR message appears the oven will be shut down.

Position of the Column Oven 4050 in Smartline Systems

Generally spoken no fixed regulations exist for positioning the Smartline Column Oven 4050 in the Smartline system. Depending on the total configuration some recommendations can be given.

To keep the dead volume as low as possible and use the shortest possible capillary the oven should be placed between the pump and the detector.

If a Smartline Autosampler 3800 or 3900 is integrated into the system the Column Oven has to be placed more on the top of the Smartline tower. Otherwise the column compartment is not sufficient accessible. In case of the Smartline Autosampler 3950 this restriction is not valid.

Installation of the Columns

Please Note the Smartline Column Oven 4050 can reach a maximum temperature of 145°C. Normal silica gel containing HPLC columns are not allowed to be operated above 90°C. Please set the allowed maximum temperature for your columns in use.

The column and capillary carrier can easily be removed from the oven compartment. Thus the capillaries and columns can be fixed outside. The whole equipped carrier is then fixed in the oven with only one screw.

Access to the Column Compartment

SOP 1 Opening and closing the column compartment

1. Open the front door at least as far as the side panel will be totally accessible.
2. Open the mechanical interlocks by screwing both black knobs counter clock wise until dead stop.
3. The side panel now can be opened and removed.
4. Before closing the compartment again make sure that both interlocks are totally opened (compare Fig. 9 A). Otherwise screw both black knobs counter clock wise until dead stop.

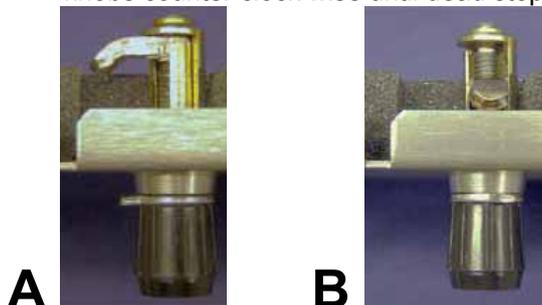
Fig. 9 Column compartment interlocks, A opened, B closed

5. Insert the side panel into the guide slots and close it.
6. Slightly pressing the side panel against the oven screw both black knobs clock wise until dead stop (Fig. 9 B).

Column and Capillary Holding

SOP 2 Installation of the capillaries

The capillaries have to be placed in the outside nuts of that holding part which is later orientated to the inner of the oven. There are four guide nuts for 1/16" OD capillaries and two for 0.5 mm OD ones (Fig. 10).

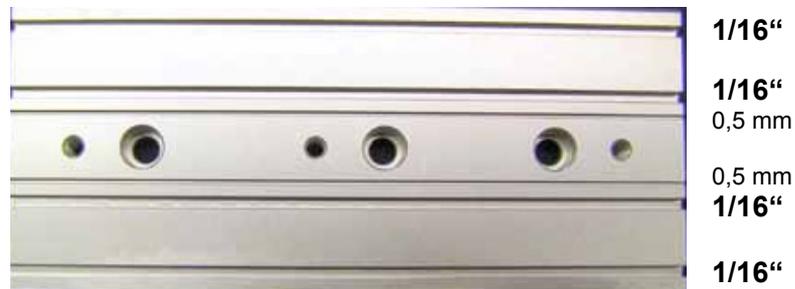
Fig. 10 Capillary guide nuts

The capillaries will be held in position with the fixing plate (Fig. 11).

Fig. 11 Fixed capillaries with 10mm screw

The screw heads of the holder serve to position the complete holder exactly during the installation of the oven (Fig. 12). They are inserted in the holes of the positioning rail.

Fig. 12 Fixed capillaries

SOP 3 Installation of the column

Insert a short screw into the middle drilling of one of the holding parts for positioning both parts to each other (Fig. 13). This is especially of advantage, if the above mentioned prefixing is not used.

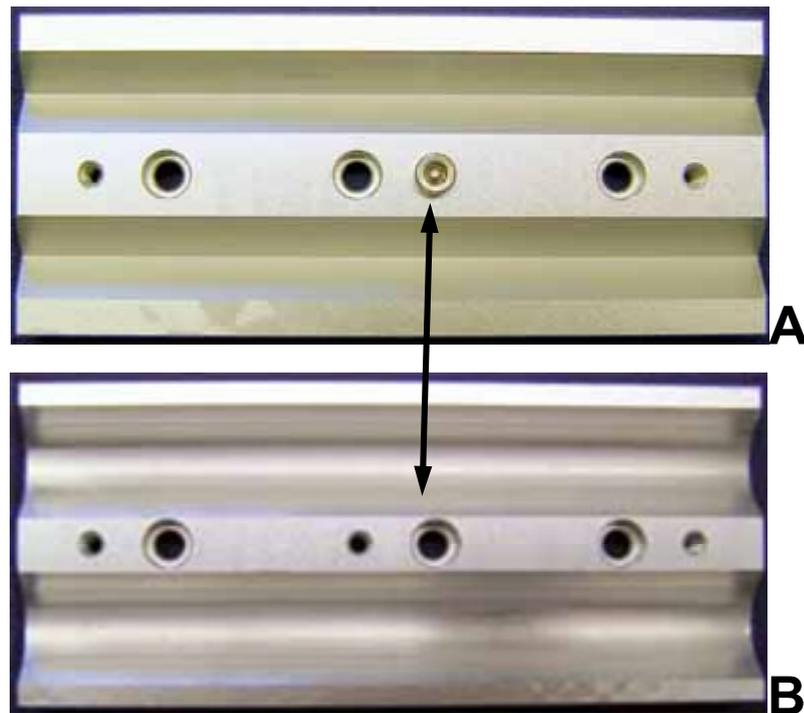
Fig. 13 Fixing the columns

For columns of different OD two types of holdings are delivered. The outer sides with their guide nuts and the drillings are for both types identically. The difference consists in the shape of the column guide to ensure as well a good fixing as an optimized heat transfer (Fig. 13 A and B).

Columns of larger diameter can be mounted either with two holding parts of type B or with one part type A and one part type B. In this case it is recommended to use the part of type A for fixing the capillaries.

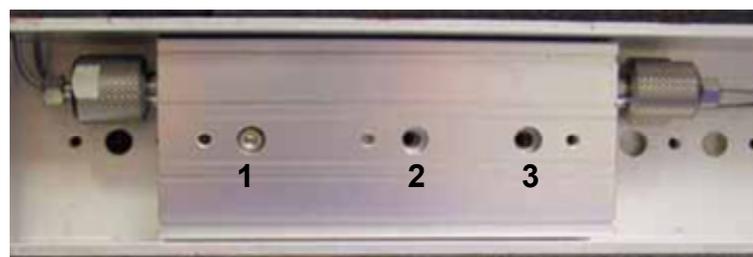
Fig. 14 Fixing the column holding in the oven

The column holder is fixated with two screws (20 mm) in the oven (with smaller holders by one screw). For the fixing use the outer drillings in the holding. (Pos. 1 and 3 in Fig. 14).

The middle drilling (Pos. 2 in Fig. 14) can be used for prefixing both the holding parts. This makes the installation of the equipped holding easier.

SOP 4 Mounting of equipped column holdings

1. Insert the holding equipped with column(s) and capillaries that way into the oven that the capillaries leave the compartment through the compartment slot and the screw heads of the capillary holding fit into two of the wholes inside the oven (Fig. 11 and Fig. 12).
2. Fix the whole holding with one screw using one of the outer drillings.

SOP 5 Connecting Capillaries

The following SOP is valid for the Smartline Column Oven 4050.

1. Push the bushing, the clamping ring, and the sealing ring onto the capillary (for standard connections just use the respective bushing and clamping ring). Please take care to note the sequence and orientation of the fittings according to Fig. 15.
2. Push the capillary as far as possible into the inlet.
3. Tighten the PEEK fitting and the Dynaseal screw fitting hand-tight. For screw fittings made of UNF metal, use a wrench to tighten.

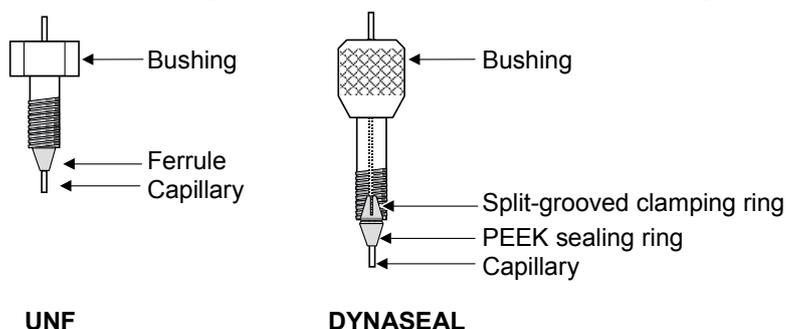
Fig. 15 Principle of the UNF- and DYNASEAL connections

To perform a post column cooling of the eluat the capillary from the column outlet is to be guided through the cooler behind the front door of the oven (Fig. 3). The cooler shall be fitted with a isolation cap to prevent humidity condensation at the cooled elements Fig. 4.

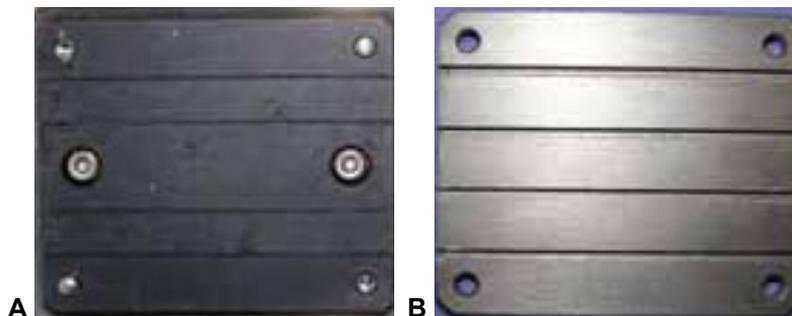
Fig. 16 Post column cooler

SOP 6 Mounting the outlet capillary into the cooler

1. Remove the isolation cap.
2. Loosen both screws at the cooler and remove the front plate.
3. Guide the capillary through one of the nuts. For 0.5 mm capillaries use a nut in the inbuilt part (Fig. 16 A) and for 1/16" capillaries a nut in the front plate (Fig. 16 B).
4. Put on the front plate and fix both screws.
5. Attach the isolation cap. Use the cuts at both sides of the cap for the capillaries.

Direct Control of the Column Oven 4050

(stand-alone-operation)

The control is performed from the main menu via the foil pad. For a not programmed thermostating only the set point must be entered and the heating switched on. To load a certain program or link enter the corresponding number and move the cursor away from the input field.

A	P--- ↑ 025°C 000 OFF [030°C]	B	L--- ↑ 025°C 000 OFF [030°C]
---	--	---	--

If the cursor is on the P (L) field you can go through and back the screens A and B using the ▲ or ▼ keys.

Input field: OFF / ON / RUN / HOLD / WAIT / REMOTE (☺)

The actual operating status of the oven is displayed on the input field in the lower left hand corner. If a program or link is loaded only the inputs „1“ and „0“ are allowed. They will start („1“ → RUN) or stop („0“ → HOLD) the program.

P05 ↑ 025°C 000 RUN 001.23 [030°C]	P05 ↑ 025°C 000 HOLD 001.74 [030°C]
--	---

Entering „0“ twice will abort the run. The oven turns to the status ON.

While a link is running also the status WAIT can be shown. The time is shown as 000.00.

L02 ↑ 040°C 000 WAIT 000.00 [045°C]

Software Control of the Column Oven 4050

Operating the Smartline Column Oven 4050 within an HPLC system controlled by one of the KNAUER software packages ChromGate® or ClarityChrom® is very easy. The Smartline Column Oven 4050 is supported in ChromGate® version 3.1.7 and higher and in all ClarityChrom® versions.

In this chapter the connections to the Column Oven 4050 necessary for working under ChromGate® or ClarityChrom® will be explained. Regarding software details refer to the corresponding instrument control manuals of the software.

If the Column Oven 4050 is under remote control the status is displayed as ☺ ON or ☹ OFF. All other displayed information corresponds to that of the stand alone operation.

P--- ↑ 025°C 000 ☹ OFF [030°C]
--

All possibilities for entering data are disabled.

Fig. 17 Remote control interfaces: LAN, RS-232 and analog connections

RS-232 Serial Interface

The RS 232 interface at the rear of the instrument (Fig. 17) enables the digital data transfer between the Column Oven 4050 and a PC with HPLC software (ChromGate® or ClarityChrom®). Connect this interface directly or occasionally using a COM port extension kit to the COM port of your PC.

LAN – Ethernet Interface

The Smartline oven S4050 can be integrated via a LAN interface into a local network. The oven gets its own IP address thereby from a DHCP server. Create your local area network (LAN) as illustrated in Fig. 18. Connect the LAN ports of the PC and all devices to be controlled using ethernet patch cables.

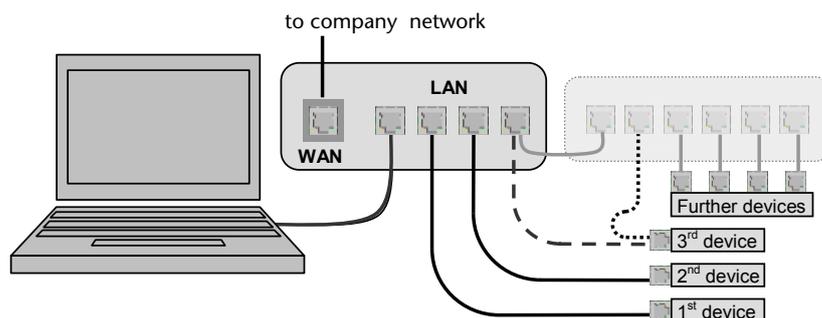
Fig. 18 Structure of the local area network (LAN)

Never use the WAN port for connection to the PC or any device to be controlled. The WAN port is exclusively to use for connection to a company network.

Make sure that only one DHCP server (normally your router) is present in your local network. All devices to be controlled must be directly connected and operate as DHCP clients (Setting: Obtain an IP address automatically). Any second DHCP server within the same network may lead to inconsistencies in IP addresses, preventing the functionality of the whole network.

To check the DHCP client setting, select *Start – Control Panel – Networking*. Right mouse click on the Ethernet icon (the annotation may differ) and then select the menu item *Properties*.

Activate in the tab *General* (possibly after scrolling) the *Internet Protocol (TCP/IP)* and click now on the *Properties* button. Compare the opening window with Fig. 19 and change the settings as far as necessary.

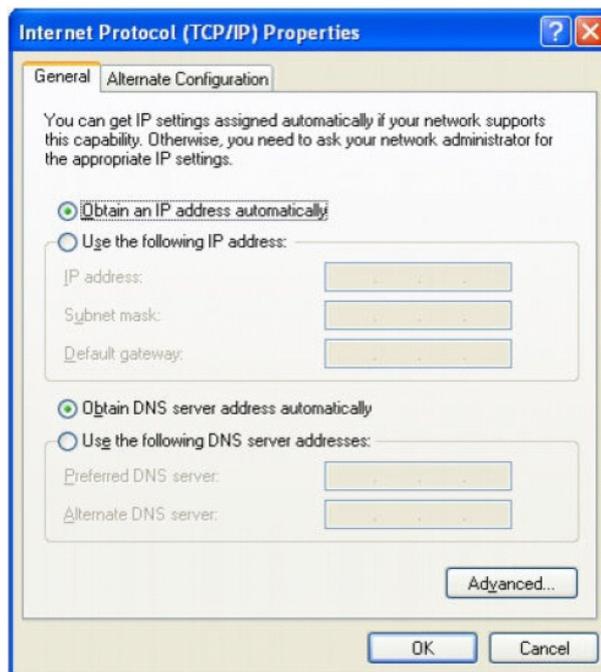
Fig. 19 Properties of the Internet Protocol (TCP/IP)

Analog Interface

Additionally the start output of an injection valve or autosampler can be connected to the START input (and GROUND) of the remote control strip.

Using the remote control socket

The remote control socket, which serves to send or receive signals from other instruments, is located on the rear panel of the Smartline Column Oven 4050 (see Fig. 20).

Please avoid touching the electrical contacts of the socket lines. Electrostatic discharges when touching the contacts could damage the electronics of the device.

Connections to the remote control socket

Two of the eight positions on the remote control socket are **ground** connections, two for ERROR communication and four serve as control connections (Start, Event):

Fig. 20 Analog remote control connectors

Control signals (events) are electric signals or pulses, which are sent during a running program via the remote control strip to control further instruments. Contacts can be closed or opened as well as 500 ms pulse signals can be activated. The event connectors are to use together with the lower GROUND connector.

Start IN	A short cut to GND switches the oven ON.
ERROR IN	An error signal remains active as long as an other instrument is in an error status.
ERROR OUT	An error signal remains active as long as the error status of the oven is indicated.
GND	This ground connector is only for Start IN, ERROR IN and ERROR OUT
EVENT 1	} Outputs to control other instruments
EVENT 2	
EVENT 3	
GND	This ground connector is only for events

Configuring Plug Strips

For external control of the Smartline Column Oven 4050 a special plug strip should be configured according to the local conditions. You can use the plug strip with 8 connectors (enclosed in the accessories) for easy installation of electrical connections.

SOP 7 Mounting Plug Strips

1. Place the connector strip (C) on a suitable surface.
2. Press down the cable terminal with the stylus (B).
3. Keep the stylus pressed down and stripped insert the cable ends (A) into the front of the connector strip.
4. Remove stylus.
5. Check whether the cables are tightly attached. Use this SOP for connecting wires to the plug strips.

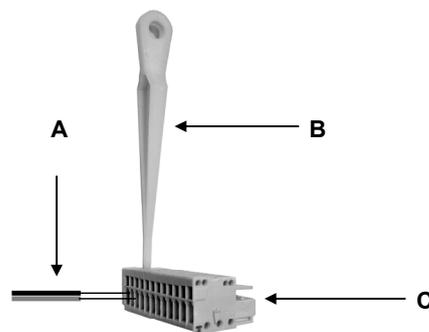
Fig. 21 Mounting of the plug strips

Installation Qualification (IQ)

Certification on the functionality of the device. During installation of the device, an installation log (IQ document) is created upon request in coordination with the technical support of the manufacturer. This installation protocol needs to be completed in full and be signed by both parties. It serves as proof of the properly executed installation and the functionality of the device.

Operation Qualification (OQ)

Extensive test of the column oven's functionality. The successfully executed OQ ensures that the column oven functions properly.

Test intervals

Run the device test at the following time intervals:

- Average use of 1 to 5 days/week: device test every 6 months
- Average use of more than 5 days/week or 24 hours/day: device test every 3 months

Execution

The execution is done either by the manufacturer's technical support or by a company authorized by the manufacturer.

The OQ documentation required for executing the OQ is with costs (once) and can be ordered separately from the manufacturer.

Maintenance and care

Proper maintenance of your HPLC device will ensure successful analyses and reproducible results.

Technical support

If you have any technical questions regarding the hardware or software of the manufacturer, please use one of the contact options below:

Hotline of the technical support department:

Languages: German and English

Available by telephone: 8 am to 5 pm (CET)

Phone: +49-(0)30-809727-0

Fax: +49-(0)30-8015010

E-mail: info@knauer.net (manufacturer)

Maintenance contract

The following maintenance work on the device may only be performed by the manufacturer or a company authorized by the manufacturer and is covered by a separate maintenance contract:

- Opening the device or removing housing parts

What maintenance tasks can users perform on the device?

Users may perform the following maintenance tasks themselves:

- Replacing the columns, holdings and fixing plates

Cleaning and caring the device

Risk of electrical shock or short circuit if cleaning solution enters the device's interior! Only lightly moisten the cleaning cloth!

All smooth surfaces of the device can be cleaned with a mild, commercially available cleaning solution, or with isopropanol.

The touchscreen of the device can be cleaned with isopropanol and wiped dry with a soft, lint-free cloth.

Environmental protection

Disposal

Drop the devices off at the local municipal waste facilities or send the devices back to the manufacturer where it will be disposed of properly.

Decontamination

Contamination of devices with toxic, infectious or radio-active substances poses a hazard for all persons during operation, repair, sale and disposal of a device.

Danger caused by toxic, infectious, or radioactive substances! A contaminated device must never be submitted for repairs, sold, or disposed of! Contract a specialist company to decontaminate the device or perform the decontamination yourself if you have the required expertise!

All contaminated devices must be properly decontaminated by a specialist company or the operating company before they can be recommissioned, repaired, sold, or disposed of.

All materials or fluids used for decontamination must be collected separately and disposed of properly.

Storage

Ambient storage conditions for the device

Air humidity: Below 90 % humidity (non-condensing)

Temperature range: 4-40°C; 39, 2-104°F

Technical Data

Ambient storage conditions	
Temperature range	4 - 40 °C; 39.2 - 104 °F
Air humidity	Below 90 % humidity (non-condensing)
Heating and cooling system	
	Micro processor controlled contact heating with thermostating of mobile phase at the column outlet and cooling at the column inlet. Active Peltier cooling of the column.
Column compartment	
Internal space	thermically isolated, shielding by air circulation
Column capacity	2 columns up to 12 mm OD and 300 mm length or 4 columns up to 12 mm OD and 125 mm length (Column fixing exchangable)
Capillary heating/cooling	
	for 1/16" and 0.5 mm
Dimensions	355 x 45 x 24 mm (inside, W x H x D)
Temperature control	
Temperature range	10°C (50°F) below room temperature – 140 °C (284°F)
Temperature stability	± 0.1 °C (32,18°F)
Technical parameters	
Display	LCD, 2 lines à 24 characters
Operation	Software controlled via digital interfaces or manuall via foil pad
Programming	Temperature gradient, 10 programms, 99 steps
Digital control and –output	Ethernet and RS-232
Leak sensor	Adjustable sensitivity, alarm
Column safety	programmable temperature shut down
GLP-Report	Serial number, Firmware version, operating time with heating/cooling, max. operating temperature
Supply voltage range	100 – 120 V / 200 – 240 V
Supply frequency	50 – 60 Hz
Power consumption	300 W
Dimensions	226 x 135 x 410 mm (W x H x D)
Weight	6.2 kg
GLP support	Service information

Delivery program

Accessories

Table 1 **Accessories kit 1**

Name	Order number
Spare rack for post column coil	A5304
Capillary holding for oven track	
Capillary 1/16 ID 0,1 x 3000	
Capillary 1/16 ID 0,1 x 3000	

Table 2 **Accessories kit 2**

Name	Order number
Spare rack for columns up to 8 mm OD	A5305
2 x half shell holdings for columns up to 8 mm OD 100 mm	
Fixing plate 100mm	
8 x ISO 4762 M3 x 10 A2 cylinder screw	
4 x ISO 4762 M3 x 20 A2 cylinder screw	
Hexagonal key no. 2,5	
Capillary 1/16 ID 0,1 x 3000	
Capillary 1/16 ID 0,1 x 3000	

Table 3 **Accessories kit 3**

Name	Order number
Spare rack for columns up to 12 mm OD	A5306
2 x half shell holdings for columns of 9 to 12 mm OD 100 mm	
Fixing plate 100mm	
8 x ISO 4762 M3 x 10 A2 cylinder screw	
4 x ISO 4762 M3 x 20 A2 cylinder screw	
Hexagonal key no. 2,5	
Capillary 1/16 ID 0,1 x 3000	
Capillary 1/16 ID 0,1 x 3000	

Table 4 **Accessories kit 4**

Name	Order number
Spare capillary set for short columns up to 125 mm length	A5307
Capillary 1/16 ID 0,1 x 3000	
Capillary 1/16 ID 0,1 x 3000	
4 x bushing UNF 10/32	
4 x ferrule 1/16" 1,4401	

Table 5 **Accessories kit 5**

Name	Order number
Spare capillary set for long columns up to 300 mm length	A5308
Capillary 1/16 ID 0,1 x 3000	
Capillary 1/16 ID 0,1 x 3000	
4 x bushing UNF 10/32	
4 x ferrule 1/16" 1,4401	

Spare parts

Table 6 **Spare parts**

Name	Order number
Power Cable type VIIG C5	M1479
Null Modem cable (complete)	A0895
10-polar flat band cable	A1467
Set WAGO Micro connect. strips for Smartline 8pol.	A1420V8
Half shell holdings for columns of 9 to 12 mm OD 40 mm	P3397
Half shell holdings for columns up to 8 mm OD 40 mm	P3399
Fixing plate 200 mm	P3447
Fixing plate 40 mm	P3401
Hexagonal key no. 2,5	X0010
Hexagonal key No. 3	X0011
Double open-end wrenches, 1/4" and 5/16"	X0003
Patch cable CAT 5e	A5255

Legal information

Warranty conditions

The factory warranty for the device is valid for 12 months after the date of dispatch. All warranty claims shall expire in the event that any unauthorized changes are made to the device.

During the warranty period, any components with material or design-related defects will be replaced or repaired by the manufacturer free of charge.

This warranty excludes the following:

- Accidental or willful damage
- Damage or errors caused by third parties that are not contractually related to the manufacturer at the time the damage occurs
- Wear parts, fuses, glass parts, columns, light sources, cuvettes and other optical components
- Damage caused by negligence or improper operation of the device and damage caused by clogged capillaries
- Packaging and transport damage

In the event of device malfunctions, contact:

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The packaging of our devices provides the best possible protection against transport damage. However, immediately inspect each delivery for signs of transport damage. If the shipment is incomplete or damaged, inform the manufacturer within three workdays. Also inform the freight carrier about transport damage.

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Smartline Column Oven 4050 order number:

A5300

complies with the following requirements and product specifications:

- DIN EN 60799 (June 1999) Electrical accessories – Cord sets and interconnection cord sets
- DIN EN 61010-1 (August 2002) Safety requirements for electrical equipment for measurement, control and laboratory use
 - Low voltage directive (2006/95/EC)
- DIN EN 61000-3-2 (March 2010) Electromagnetic compatibility (EMC) Part 3-2
EN 61000-3-2:2006 + A1:2009 + A2:2009
 - EMC directive (2004/108/EC)
- DIN EN 61326-1 (October 2006) Electrical equipment for measurement, control and laboratory use – EMC requirements
 - DIN EN 61326-1 Corrigendum 2 (April 2011)
- Directives for an environmentally sound use of electrical and electronic equipment
 - RoHS directive 2002/95/EC (February 2003) on the restriction of the use of certain hazardous substances in electrical and electronic equipment
 - WEEE directive 2002/96/EC (February 2003) on waste electrical and electronic equipment

Berlin, 2011-08-04

Dr. Alexander Bünz (Managing Director)

The mark of conformity has been applied to the rear panel of the device.

Index

- A**
- Ambient conditions 32, 33
 - Analog interface..... 29
 - Arrow keys..... 14
 - Auto test..... 16
 - Autosampler 23
- C**
- Care 31
 - ChromGate 27, 48, 63
 - ClarityChrom..... 27, 63
 - Cleaning..... 32
 - Column 23
- D**
- Declaration of conformity..... 37
 - Decontamination..... 32
 - Device test..... 31
 - Installation Qualification..... 31
 - Device types 6
 - DHCP-Client 28
 - Direct control 27
 - Disposal..... 32
 - DYNASEAL 26
- E**
- Editing menu..... 19
 - END - key 14
 - Environmental protection..... 32
 - Error screen 22
 - Ethernet 28
 - Events field 17
 - External events..... 17
- F**
- Flammability 7
 - Foil keys..... 14
 - Front elements..... 15
 - Front view 13
 - Function elements 13
- G**
- General description 12
 - GLP menu 22
- H**
- Hysteresis..... 19
- I**
- Installation 10
 - Installation site..... 11
 - Intended Use 6
 - Device types..... 6
 - Interface..... 18
 - Internal software 17
- L**
- Laboratory use 6
 - LAN 28
 - Leak sensor..... 18
 - Lieferprogramm..... 34
 - Link..... 21
- M**
- Mains connection 8
 - Maintenance..... 31
 - Maintenance contract..... 31
 - Manufacturer 36
- O**
- ON OFF - key..... 14
 - On/off switch 17
 - Operating voltage..... 16
 - Operation
 - Device 8
 - Original accessories..... 10
- P**
- PEEK..... 8
 - Plug strips 30
 - Post column cooler 15, 26
 - Power supply..... 8, 16
 - Program 19
 - Protective measures 8
- R**
- Rear pannel elements 16
 - Rearview 16
 - Remote control..... 29
 - RS-232 interface 28
 - Run time 20
- S**
- Safety 6
 - Scope of delivery 10
 - Self-ignition point 7
 - Set point 17
 - Setup..... 18
 - Smartline System..... 23
 - Software Control 27
 - Software structure..... 17
 - Solvent
 - suitable..... 7
 - toxicity 7
 - unsuitable..... 7
 - Solvent tray 7
 - Solvents 7
 - SOP's..... 5

Space requirements	11
Spare parts	35
Stand alone	27
Standby key	14
Storage	32
Symbols	9
System test	17

T

Target group	8
Technical Data	33
Technical support	31
Temperature	18

Toxicity	
Solvent	7

U

Unpacking	13
-----------------	----

W

Wait status	21
WAN	28
Warnings	9
Warranty conditions	36

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All rights reserved.
Dr. Ing. Herbert Knauer GmbH
Dr. Ing. Herbert Knauer GmbH Telefax: +49-(0)30-8015010